

## Folders should be organized:


- 1. Overarching folder (ex. 3HAA potatoes)
  - a. Repeats (ex. Repeat\_1)
    - i. Full condition name (ex. GLS130 control)
      - 1. Days (...GLS130 control D1)
- 1. Run SICKO\_2022
- 2. Select Full\_Condition\_Name
- 3. WARNING box will prompt
  - a. If continue selected, all csv files in the folder will be deleted. Save data elsewhere if you don't want them to be deleted.

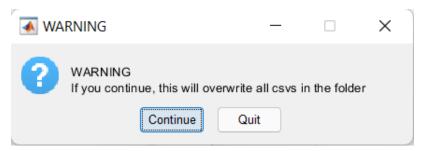


b. If quit selected, program will quit

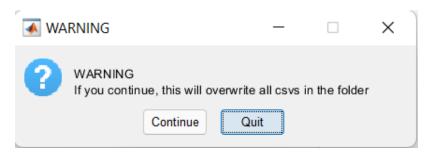

## 4. GUI will prompt

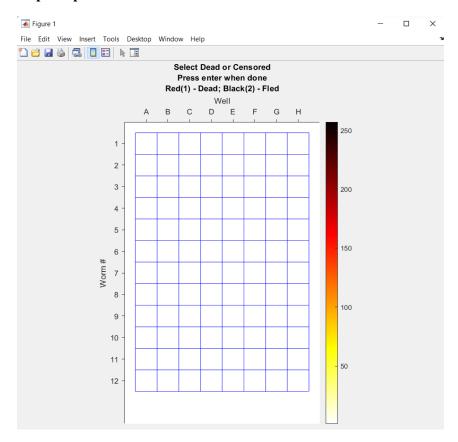


a. Select worms that are dead (Red, 1 click) and fled (Black, 2 clicks)

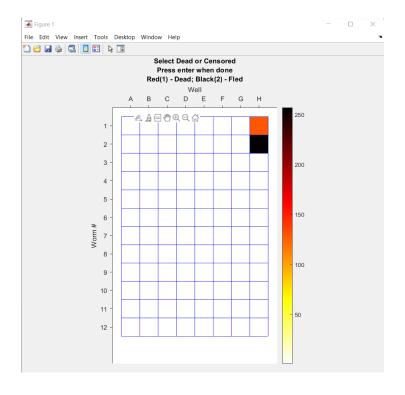


- b. Double click or click enter to confirm
- c. Save box will prompt

## d. Select Y to save

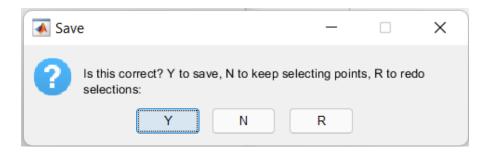


e. Select N to keep selecting points

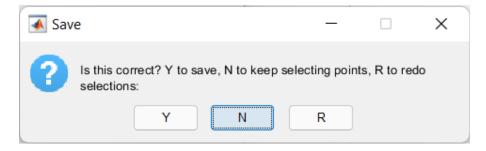


f. Select R to erase previous selections and select again




5. Enter the day that the worm died or fled

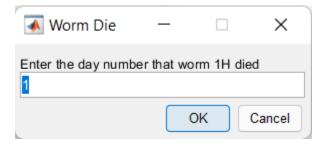


- 6. Beginning Data Processing
- 7. Redo? With prompt

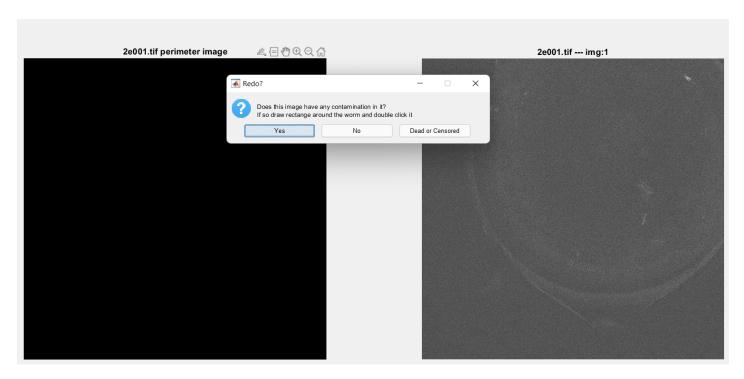


a. Select Yes to draw a rectangle around the worm and double click to confirm.

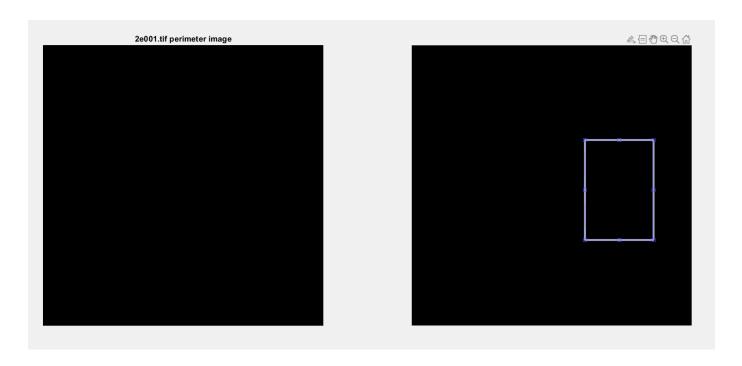


b. Select No to continue to the next photo

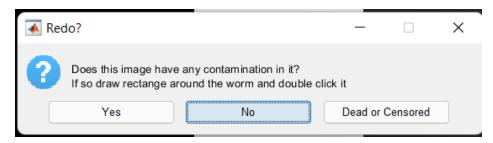


c. Select Dead or Censored for more options

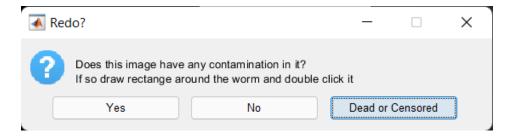


d. Select Dead to declare worm dead, crop for contamination will prompt

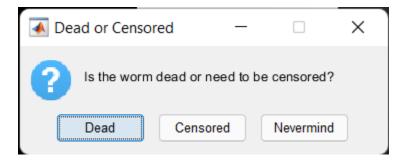

e. Select Yes to crop for contamination, Select No to move to the next image

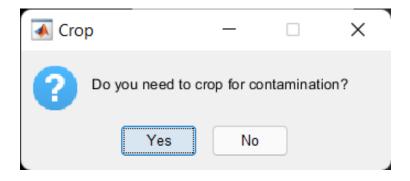

f. If Yes, draw a rectangle around the portion of the image you want to keep and double click

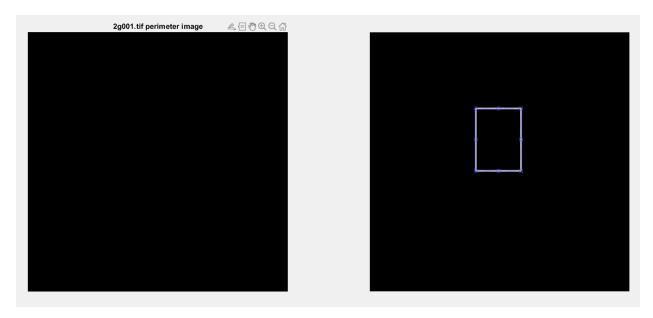

g. Select Censored to censor the worm if image cannot be used



h. Select Nevermind to move to the next





- if you want to exit the program
- 9. Select Yes to quit, Select No to continue

